

Instrument controller software run summary:

Filename and data path: C:\Agilent Technologies\Data\2021 06 21\13-15-04\2021 06 21 13H 15M.raw

Created: Monday, June 21, 2021 1:36:48 PM

Number of capillaries: 3

Array serial number: 091520-10SFS

Effect length: 33 cm

Array usage count: 329

Instrument type: 5300 Fragment Analyzer

Instrument controller software version: 3.1.0.12

Device serial number: 3178

Method Information

Method name: DNF-474-33 - HS NGS Fragment 1-6000bp.mthds

Gel prime: No

Full conditioning: Yes

Gel prime to buffer: No

Gel selection: Gel 1

Perform prerun: 6.0 kV, 30 sec.

Rinse: No

Marker 1: No

Rinse: Tray: 3, Row: A, Dip count: 1

Sample injection: 5.0 kV, 30 sec.

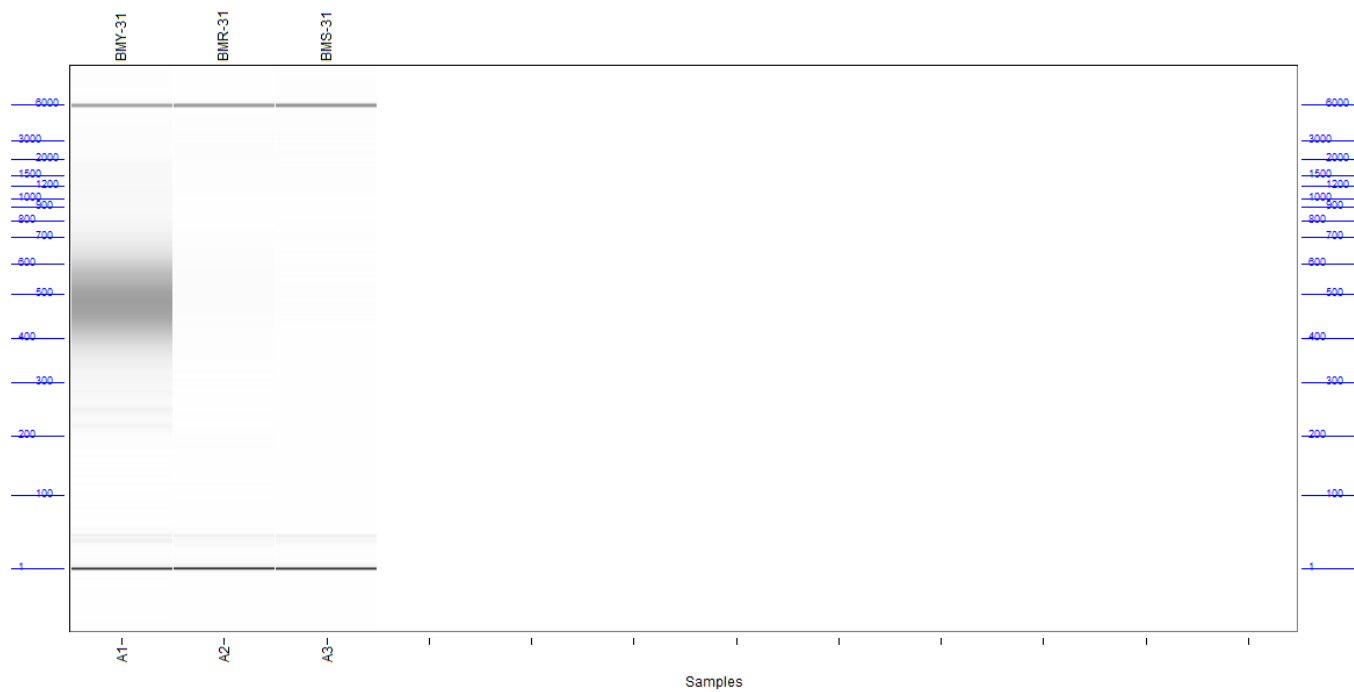
Separation: 6.0 kV, 50.0 min.

Tray name: Tray-1

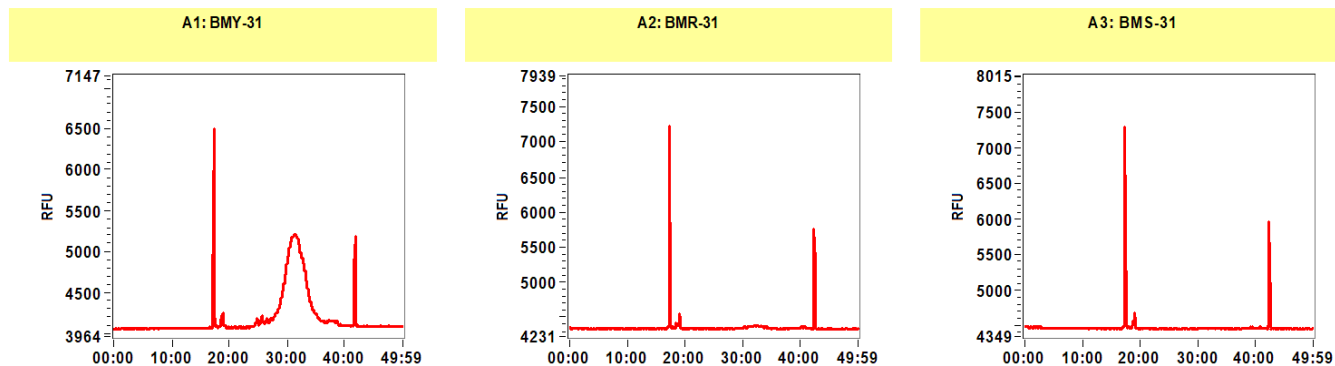
Analysis mode: NGS

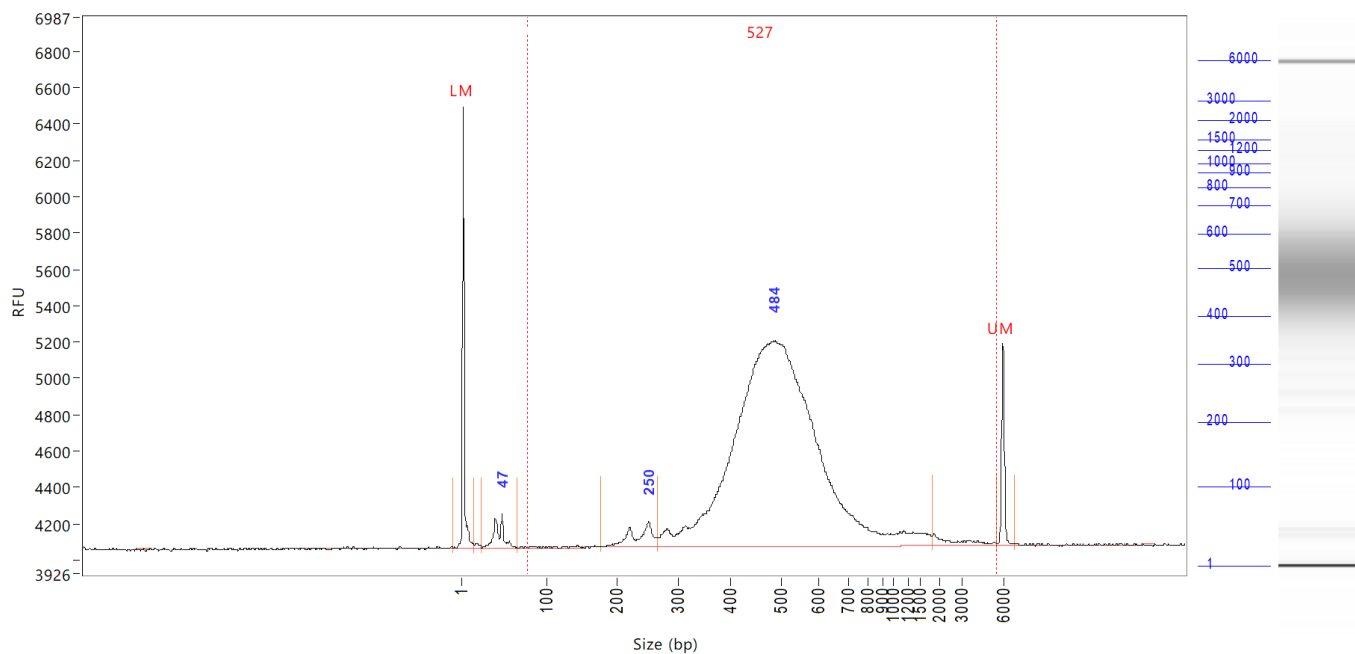
Notes

Gel Image



Filename and data path: C:\Agilent Technologies\Data\2021 06 21\13-15-04\2021 06 21 13H 15M.raw



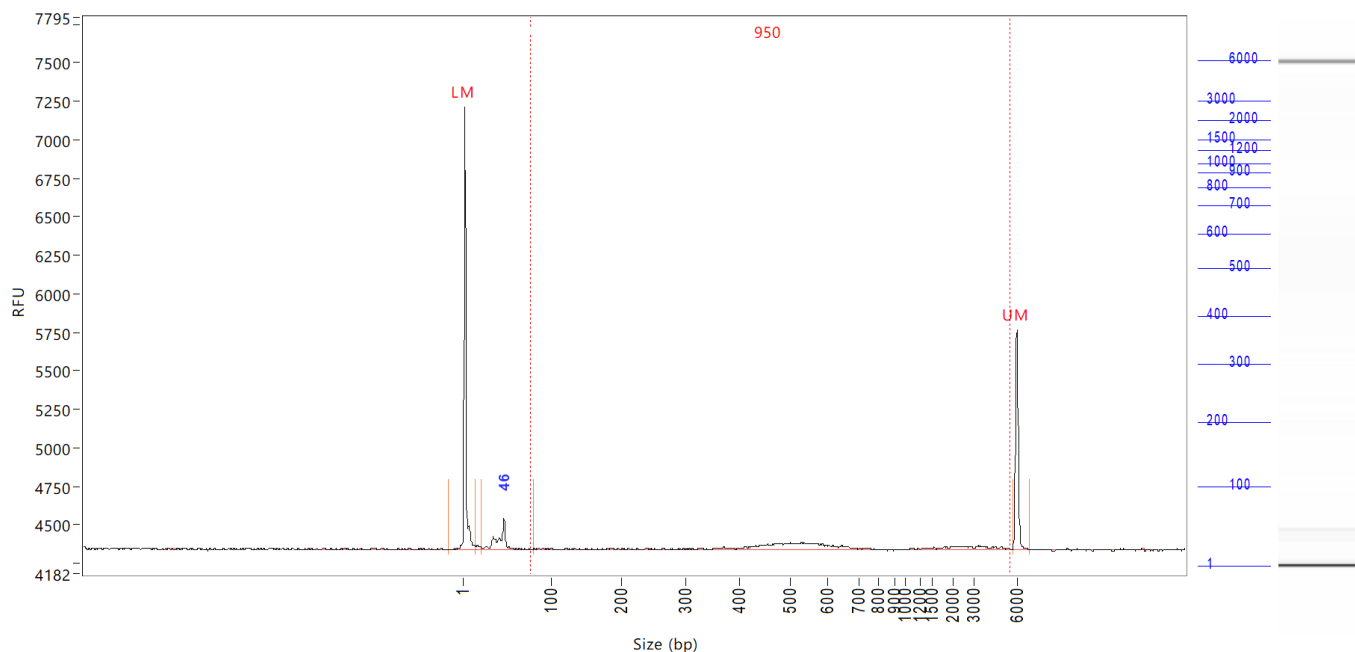
Sample: BMY-31**Well location:** A1**Created:** Monday, June 21, 2021 1:36:48 PM

Peak	Size	Concentration	From	To	Average size	CV%	RFU	Corrected peak area
	(bp)	(ng/uL)	(bp)	(bp)	(bp)			
1	1 (LM)	0.0174	0	15	1	202.02	2431	13.671
2	47	0.0613	23	65	42	15.66	187	4.005
3	250	0.0777	178	264	233	8.59	141	5.078
4	484	2.4095	264	1812	511	31.29	1132	157.409
5	6000 (UM)	0.0050	5553	6784	5990	2.10	1115	3.936

TIC: 2.5486 ng/uL
 TIM: 10.6961 nmole/L
 Total concentration: 2.5879 ng/uL

Smear Analysis 75 bp to 5500 bp 2.5210 ng/uL 97.4 %Total 7.8665 nmole/L 527 Avg. Size (bp) 61.02 %CV

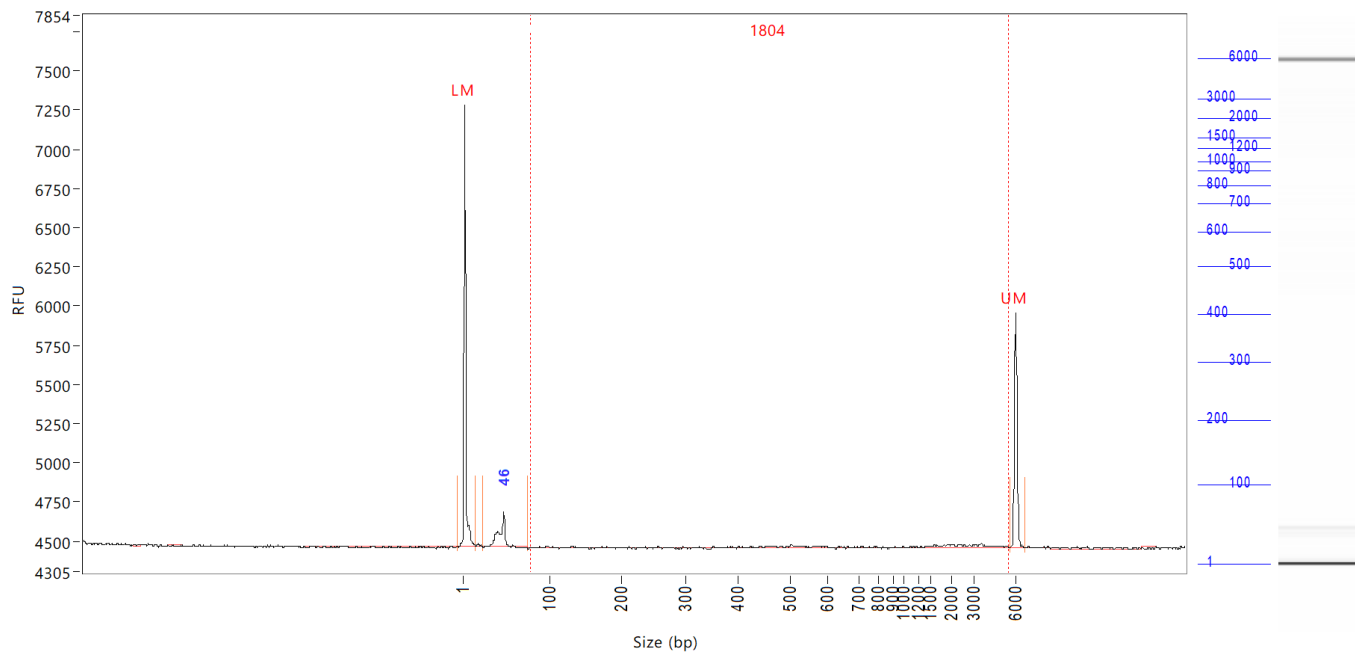
Sample peak width (sec): 50 Sample min peak height: 25 Sample baseline V to V?: Y Sample baseline V to V points: 3
 Sample filter: Binomial Number of points for filter: 3 Sample start region (min): 0 Sample end region (min): 50
 Manual baseline start (min): 10 Manual baseline end (min): 48
 Marker peak width (sec): 5 Marker min peak height: 200 Marker baseline V to V?: Y Marker baseline V to V points: 3
 Lower marker selection: First peak > 200 RFU Upper marker selection: Last peak > 200 RFU
 Ladder size (bp): 1, 100, 200, 300, 400, 500, 600, 700, 800, 900, 1000, 1200, 1500, 2000, 3000, 6000
 Quantification using: Ladder Final concentration (ng/uL): 0.0830 Dilution factor: 12.0

Sample: BMR-31**Well location:** A2**Created:** Monday, June 21, 2021 1:36:48 PM

Peak	Size	Concentration	From	To	Average size	CV%	RFU	Corrected peak area
	(bp)	(ng/uL)	(bp)	(bp)	(bp)			
1	1 (LM)	0.0174	0	15	1	177.43	2877	16.086
2	46	0.0449	22	81	41	20.56	206	3.449
3	6000 (UM)	0.0054	5684	6941	5995	2.08	1426	4.973
	TIC:	0.0449	ng/uL					
	TIM:	1.7906	nmole/L					
	Total concentration:	0.1560	ng/uL					

Smear Analysis 75 bp to 5500 bp 0.1074 ng/ul 68.9 %Total 0.1860 nmole/L 950 Avg. Size (bp) 113.61 %CV

Sample peak width (sec): 50 Sample min peak height: 25 Sample baseline V to V?: Y Sample baseline V to V points: 3
Sample filter: Binomial Number of points for filter: 3 Sample start region (min): 0 Sample end region (min): 50
Manual baseline start (min): 10 Manual baseline end (min): 48
Marker peak width (sec): 5 Marker min peak height: 200 Marker baseline V to V?: Y Marker baseline V to V points: 3
Lower marker selection: First peak > 200 RFU Upper marker selection: Last peak > 200 RFU
Ladder size (bp): 1, 100, 200, 300, 400, 500, 600, 700, 800, 900, 1000, 1200, 1500, 2000, 3000, 6000
Quantification using: Ladder Final concentration (ng/uL): 0.0830 Dilution factor: 12.0

Sample: BMS-31**Well location:** A3**Created:** Monday, June 21, 2021 1:36:48 PM

Peak	Size	Concentration	From	To	Average size	CV%	RFU	Corrected peak area
	(bp)	(ng/uL)	(bp)	(bp)	(bp)			
1	1 (LM)	0.0174	0	14	1	184.62	2823	15.434
2	46	0.0448	23	73	41	12.76	223	3.308
3	6000 (UM)	0.0058	5553	6653	5991	1.80	1498	5.126
	TIC:	0.0448	ng/uL					
	TIM:	1.7895	nmole/L					
	Total	0.0860	ng/uL					
	concentration:							

Smear Analysis 75 bp to 5500 bp 0.0387 ng/ul 45.1 %Total 0.0353 nmole/L 1804 Avg. Size (bp) 78.05 %CV

Sample peak width (sec): 50 Sample min peak height: 25 Sample baseline V to V?: Y Sample baseline V to V points: 3
Sample filter: Binomial Number of points for filter: 3 Sample start region (min): 0 Sample end region (min): 50
Manual baseline start (min): 10 Manual baseline end (min): 48
Marker peak width (sec): 5 Marker min peak height: 200 Marker baseline V to V?: Y Marker baseline V to V points: 3
Lower marker selection: First peak > 200 RFU Upper marker selection: Last peak > 200 RFU
Ladder size (bp): 1, 100, 200, 300, 400, 500, 600, 700, 800, 900, 1000, 1200, 1500, 2000, 3000, 6000
Quantification using: Ladder Final concentration (ng/uL): 0.0830 Dilution factor: 12.0

Sample: ladder**Well location:** D12**Created:** Monday, June 21, 2021 1:36:48 PM**Fit type:** Point to point

Calibration curve

